



Rapid Insulin Test System Product Code: 5825-300

1.0 INTRODUCTION

Monobind Insulin Microplate Elisa test is intended to be used for the quantitative determination of insulin levels in human serum.

2.0 SUMMARY AND EXPLANATION OF THE TEST

Human insulin is a peptide produced in the beta cells of the pancreas and is responsible for the metabolism and storage of carbohydrates. As a result of biofeedback the insulin levels increase with intake of sugars and decline when sugar content is low for absorption. In the diabetic population the mechanism of insulin production is impaired because of genetic predispositions (**Type I**) or because of lifestyle and/or hereditary factors (**Type II**). In such cases either the insulin production has to be boosted by medication or it has to be supplemented by oral or intravenous methods. The quantitative determination of insulin can help in dose selection the patient has to be subjected to.


On the other hand the circulatory insulin can be found at much higher levels in patients with pancreatic tumors. These tumors secrete abnormally high levels of insulin and thus cause hypoglycemia. Accordingly, fasting hypoglycemia associated with inappropriately high concentrations of insulin strongly suggests an islet-cell tumor (insulinoma). To distinguish insulinomas from factitious hypoglycemia due to insulin administration, serum C-peptide values are recommended. (*Please see Monobind C-Peptide Microwell Elisa Cat#2525-300*). These insulinomas can be localized by provocative intravenous doses of *tolbutamide* and *calcium*.

3.0 PRINCIPLE

Immunoenzymometric assay (TYPE 3):

The essential reagents required for an immunoenzymometric assay include high affinity and specificity antibodies (Ab), (enzyme conjugated and immobilized), with different and distinct epitope recognition, **in excess**, and native antigen (Ag). In this procedure, the immobilization takes place during the assay at the surface of a microplate well through the interaction of streptavidin coated on the well and exogenously added biotinylated monoclonal Insulin antibody.

A reaction results upon mixing monoclonal biotinylated antibody, the enzyme-labeled antibody and a serum containing the native antigen, without competition or steric hindrance, to form a soluble sandwich complex. The interaction is illustrated by the following equation:

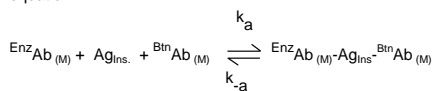

$\text{B}^{\text{In}}\text{Ab}_{(M)}$ = Biotinylated Monoclonal Ab (Excess Quantity)

Ag_{Ins} = Native Antigen (Variable Quantity)

$\text{EnzAb}_{(M)}$ = Enzyme labeled Monoclonal Ab (Excess Quantity)

$\text{EnzAb}_{(M)}\text{-Ag}_{\text{Ins}}\text{-B}^{\text{In}}\text{Ab}_{(M)}$ = Antigen-Antibodies complex

k_a = Rate Constant of Association

k_a = Rate Constant of Dissociation

Simultaneously, the complex is deposited to the well through the high affinity reaction of streptavidin and biotinylated antibody. This interaction is illustrated below:



$\text{Streptavidin}_{\text{CW}}$ = Streptavidin immobilized on well
 $\text{Immobilized complex}$ = sandwich complex bound to the solid surface

After equilibrium is attained, the antibody-bound fraction is separated from unbound antigen by decantation or aspiration. The enzyme activity in the antibody-bound fraction is directly proportional to the native antigen concentration. By utilizing several different serum references of known antigen values, a dose response curve can be generated from which the antigen concentration of an unknown can be ascertained.

4.0 REAGENTS

Materials Provided:

A. Insulin Calibrators – 2.0 ml/vial (Dried) – [Icons A – F]

Six (6) vials of references for Insulin antigen at levels of 0(A), 5(B), 25(C), 50(D), 100(E), and 300(F) $\mu\text{U/ml}$. Reconstitute each vial with 2ml of distilled or deionized water. The reconstituted calibrators are stable for sixty (60) days at 2-8°C. A preservative has been added.

Note: The calibrators, human serum based, were calibrated using a reference preparation, which was assayed against the WHO 1st IRP 66/304.

B. Insulin Enzyme Reagent – 13ml/vial - Icon

One (1) vial containing enzyme labeled affinity purified monoclonal mouse x-insulin IgG, biotinylated monoclonal mouse x-insulin IgG in buffer, dye, and preservative. Store at 2-8°C.

C. Streptavidin Coated Plate – 96 wells – Icon

One 96-well microplate coated with streptavidin and packaged in an aluminum bag with a drying agent. Store at 2-8°C.

D. Wash Solution Concentrate – 20 ml [Icon

One (1) vial contains a surfactant in buffered saline. A preservative has been added. Store at 2-8°C.

E. Substrate Solution – 12ml/vial - Icon

One (1) bottle contains tetramethylbenzidine (TMB) and hydrogen peroxide (H_2O_2) in buffer. Store at 2-8°C.

F. Stop Solution – 8ml/vial - Icon

One (1) bottle contains a strong acid (1N H_2SO_4). Store at 2-30°C.

G. Product Instructions

Note 1: Do not use reagents beyond the kit expiration date.

Note 2: Avoid extended exposure to heat and light. **Opened reagents are stable for sixty (60) days when stored at 2-8°C. Kit and component stability are identified on the label.**

Note 3: Above reagents are for a single 96-well microplate.

4.1 Required But Not Provided:

- Pipette(s) capable of delivering 50 μl and 100 μl volumes with a precision of better than 1.5%.
- Dispenser(s) for repetitive deliveries of 0.100ml and 0.350ml volumes with a precision of better than 1.5% (optional).
- Microplate washer or a squeeze bottle (optional).
- Microplate Reader with 450nm and 620nm wavelength absorbance capability (The 620nm filter is optional)
- Absorbent Paper for blotting the microplate wells.
- Plastic wrap or microplate cover for incubation steps.
- Vacuum aspirator (optional) for wash steps.
- Timer.
- Storage container for storage of wash buffer.
- Distilled or deionized water.
- Quality Control Materials.

5.0 PRECAUTIONS

**For In Vitro Diagnostic Use
Not for Internal or External Use in Humans or Animals**

All products that contain human serum have been found to be non-reactive for Hepatitis B Surface antigen, HIV 1&2 and HCV antibodies by FDA required tests. Since no known test can offer complete assurance that infectious agents are absent, all human serum products should be handled as potentially hazardous and capable of transmitting disease. Good laboratory procedures for handling blood products can be found in the Center for Disease Control / National Institute of Health, "Biosafety in Microbiological and Biomedical Laboratories," 2nd Edition, 1988, HHS.

Safe Disposal of kit components must be according to local regulatory and statutory requirement.

6.0 SPECIMEN COLLECTION AND PREPARATION

The specimens shall be blood; serum in type, and the usual precautions in the collection of venipuncture samples should be observed. For accurate comparison to established normal values, a fasting morning serum sample should be obtained. The blood should be collected in a plain redtop venipuncture tube without additives or anti-coagulants (for serum) or evacuated tube(s) containing EDTA or heparin. Allow the blood to clot for serum samples. Centrifuge the specimen to separate the serum or plasma from the cells.

Samples may be refrigerated at 2-8°C for a maximum period of five (5) days. If the specimen(s) cannot be assayed within this time, the sample(s) may be stored at temperatures of -20°C for up to 30 days. Avoid use of contaminated devices. Avoid repetitive freezing and thawing. When assayed in duplicate, 0.100ml of the specimen is required.

7.0 QUALITY CONTROL

Each laboratory should assay controls at levels in the low, normal and elevated range for monitoring assay performance. These controls should be treated as unknowns and values determined in every test procedure performed. Quality control charts should be maintained to follow the performance of the supplied reagents. Pertinent statistical methods should be employed to ascertain trends. Significant deviation from established performance can indicate unnoticed change in experimental conditions or degradation of kit reagents. Fresh reagents should be used to determine the reason for the variations.

8.0 REAGENT PREPARATION:

1. Wash Buffer

Dilute contents of wash concentrate to 1000ml with distilled or deionized water in a suitable storage container. Diluted buffer can be stored at 2-30°C for up to 60 days.

Note 1: Do not use the working substrate if it looks blue.

Note 2: Do not use reagents that are contaminated or have bacteria growth.

9.0 TEST PROCEDURE

Before proceeding with the assay, bring all reagents, serum references and controls to room temperature (20 - 27°C).

****Test Procedure should be performed by a skilled individual or trained professional****

- Format the microplates' wells for calibrator, control and patient specimen to be assayed in duplicate. **Replace any unused microwell strips back into the aluminum bag, seal and store at 2-8°C.**
- Pipette 0.050 ml (50 μl) of the appropriate calibrators, controls and samples into the assigned wells.
- Add 0.100 ml (100 μl) of the Insulin Enzyme Reagent to each well. **It is very important to dispense all reagents close to the bottom of the microwell.**
- Swirl the microplate gently for 20-30 seconds to mix. Cover with a plastic wrap.
- Incubate for 60 minutes at room temperature (20-27°C).
- Discard the contents of the microplate by decantation or aspiration. If decanting, tap and blot the plate dry with absorbent paper.

- Add 350 μl of wash buffer (see Reagent Preparation Section), decant (tap and blot) or aspirate. Repeat two (2) additional times for a total of three (3) washes. **An automatic or manual plate washer can be used. Follow the manufacturer's instruction for proper usage. If a squeeze bottle is used, fill each well to the top by squeezing the container. Avoiding air bubbles. Decant the wash and repeat two (2) additional times.**
- Add 0.100 ml (100 μl) of substrate solution to all wells.

DO NOT SHAKE THE PLATE AFTER SUBSTRATE ADDITION

- Incubate at room temperature for fifteen (15) minutes.
- Add 0.050ml (50 μl) of stop solution to each well and mix gently for 15-20 seconds. **Always add reagents in the same order to minimize reaction time differences between wells.**
- Read the absorbance in each well at 450nm (using a reference wavelength of 620-630nm to minimize well imperfections) in a microplate reader. **The results should be read within fifteen (15) minutes of adding the stop solution.**

10.0 CALCULATION OF RESULTS

A dose response curve is used to ascertain the concentration of Insulin in unknown specimens.


- Record the absorbance obtained from the printout of the microplate reader as outlined in Example 1.
- Plot the absorbance for each duplicate serum reference versus the corresponding Insulin concentration in $\mu\text{U/ml}$ on linear graph paper (do not average the duplicates of the serum references before plotting).
- Draw the best-fit curve through the plotted points.
- To determine the concentration of Insulin for an unknown, locate the average absorbance of the duplicates for each unknown on the vertical axis of the graph, find the intersecting point on the curve, and read the concentration (in $\mu\text{U/ml}$) from the horizontal axis of the graph (the duplicates of the unknown may be averaged as indicated). In the following example, the average absorbance (0.514) intersects the dose response curve at 45.5 $\mu\text{U/ml}$ for the Insulin concentration (See Figure 1).

Note: Computer data reduction software designed for ELISA assays may also be used for the data reduction. **If such software is utilized, the validation of the software should be ascertained.**

EXAMPLE 1

Sample I.D.	Well Number	Abs (A)	Mean Abs (B)	Value ($\mu\text{U/ml}$)
Cal A	A1	0.019	0.019	0
	B1	0.019		
Cal B	C1	0.071	0.071	5
	D1	0.070		
Cal C	E1	0.284	0.296	25
	F1	0.311		
Cal D	G1	0.548	0.562	50
	H1	0.576		
Cal E	A2	1.155	1.160	100
	B2	1.166		
Cal F	C2	2.635	2.694	300
	D2	2.752		
Unknown	E2	0.513	0.514	45.5
	F2	0.515		

*The data presented in Example 1 and Figure 1 are for illustration only and **should not** be used in lieu of a dose response curve prepared with each assay.

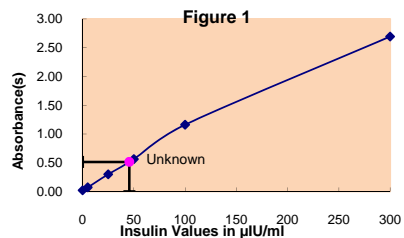

11.0 QC PARAMETERS

In order for the assay results to be considered valid the following criteria should be met:

- The OD* of 0 µIU/ml calibrator should be ≤ 0.1 .
- The OD of the 300µIU/ml calibrator should be ≥ 1.3 .
- Four out of six quality control pools should be within the established ranges.

*OD = Optical Density (Absorbance at 450nm)

12.0 RISK ANALYSIS

The MSDS and Risk Analysis Form for this product is available on request from Monobind Inc.

12.1 Assay Performance

- It is important that the time of reaction in each well is held constant to achieve reproducible results.
- Pipetting of samples should not extend beyond ten (10) minutes to avoid assay drift.
- Highly lipemic, hemolyzed or grossly contaminated specimen(s) should not be used.
- If more than one (1) plate is used, it is recommended to repeat the dose response curve.
- The addition of substrate solution initiates a kinetic reaction, which is terminated by the addition of the stop solution. Therefore, the substrate and stop solution should be added in the same sequence to eliminate any time-deviation during reaction.
- Plate readers measure vertically. Do not touch the bottom of the wells.
- Failure to remove adhering solution adequately in the aspiration or decantation wash step(s) may result in poor replication and spurious results.
- Use components from the same lot. No intermixing of reagents from different batches.
- Patient samples with Insulin concentrations above 300µIU/ml may be diluted with the zero calibrator and re-assayed. Multiply the value obtained by the dilution factor to obtain the corrected value.
- Accurate and precise pipetting, as well as following the exact time and temperature requirements prescribed are essential. Any deviation from Monobind's IFU may yield inaccurate results.
- All applicable national standards, regulations and laws, including, but not limited to, good laboratory procedures, must be strictly followed to ensure compliance and proper device usage.
- It is important to calibrate all the equipment e.g. Pipettes, Readers, Washers and/or the automated instruments used with this device, and to perform routine preventative maintenance.
- Risk Analysis- as required by CE Mark IVD Directive 98/79/EC - for this and other devices, made by Monobind, can be requested via email from Monobind@monobind.com.

12.2 Interpretation

- Measurements and interpretation of results must be performed by a skilled individual or trained professional.
- Laboratory results alone are only one aspect for determining patient care and should not be the sole basis for therapy, particularly if the results conflict with other determinants.
- For valid test results, adequate controls and other parameters must be within the listed ranges and assay requirements.
- If test kits are altered, such as by mixing parts of different kits, which could produce false test results, or if results are incorrectly interpreted, **Monobind shall have no liability.**

- If computer controlled data reduction is used to interpret the results of the test, it is imperative that the predicted values for the calibrators fall within 10% of the assigned concentrations.

13.0 EXPECTED VALUES

Insulin values are consistently higher in plasma than in serum; thus, serum is preferred. Compared with fasting values in non-obese non-diabetic individuals, insulin levels are higher in obese non-diabetic subjects and lower in trained athletes. Although proinsulin cross reacts with most competitive insulin assays, there is less than 1% cross reaction found with proinsulin using Monobind Insulin Microwell Elisa.

Based on the clinical data gathered by Monobind in concordance with the published literature the following ranges have been assigned. **These ranges should be used as guidelines only:**

POPULATION	RANGE
Children < 12 yrs	< 10 µIU/ml
Adult (Normal)	0.7 – 9.0 µIU/ml
Diabetic (Type II)	0.7 – 25 µIU/ml

It is important to keep in mind that establishment of a range of values which can be expected to be found by a given method for a population of "normal"-persons is dependent upon a multiplicity of factors: the specificity of the method, the population tested and the precision of the method in the hands of the analyst. For these reasons each laboratory should depend upon the range of expected values established by the Manufacturer only until an in-house range can be determined by the analysts using the method with a population indigenous to the area in which the laboratory is located.

14.0 PERFORMANCE CHARACTERISTICS

14.1 Precision

The within and between assay precision of the Insulin AccuBind™ ELISA Microplate Test System were determined by analyses on three different levels of pool control sera. The number, mean value, standard deviation and coefficient of variation for each of these control sera are presented in Table 2 and Table 3.

TABLE 2
Within Assay Precision (Values in µIU/ml)

Sample	N	X	σ	C.V.
Pool 1	20	11.6	0.93	8.0%
Pool 2	20	39.8	1.96	4.9%
Pool 3	20	117.8	6.00	5.1%

TABLE 3
Between Assay Precision (Values in µIU/ml)

Sample	N	X	σ	C.V.
Pool 1	10	10.3	1.01	9.8%
Pool 2	10	41.5	2.32	5.6%
Pool 3	10	122.1	8.32	6.8%

*As measured in ten experiments in duplicate over seven days.

14.2 Sensitivity

The sensitivity (detection limit) was ascertained by determining the variability of the 0 µIU/ml serum calibrator and using the 2σ (95% certainty) statistic to calculate the minimum dose. The assay sensitivity was found to be 0.75 µIU/ml.

14.3 Accuracy

The Monobind Insulin AccuBind™ ELISA was compared with a reference coated tube radioimmunoassay assay. Biological specimens from population (symptomatic and asymptomatic) were used. (The values ranged from 0.01µIU/ml – 129µIU/ml). The total number of such specimens was 104. The data obtained is displayed in Table 4.

TABLE 4

Method	Mean (x)	Least Square Regression Analysis	Correlation Coefficient
This Method	13.6	$y = 2.6 + 0.91(x)$	0.975
Reference	11.4		

Only slight amounts of bias between the Insulin AccuBind™ ELISA system and the reference method are indicated by the closeness of the mean values. The least square regression

equation and correlation coefficient indicates excellent method agreement.

14.4 Specificity

The cross-reactivity of the Insulin AccuBind™ ELISA method to selected substances was evaluated by adding the interfering substance(s) to a serum matrix at the following concentration(s). The cross-reactivity was calculated by deriving a ratio between dose of interfering substance to dose of insulin needed to produce the same absorbance.

Substance	Cross Reactivity	Concentration
Insulin	1.0000	-
Proinsulin	0.0078	100 ng/ml
C-Peptide	N/D	75 ng/ml
Glucagon	N/D	150 ng/ml

15.0 REFERENCES

- Eastham RD, *Biochemical Values in Clinical Medicine*, 7th Ed Bristol, England, John Wright & Sons, Ltd (1985).
- Gerbitz VKD, "Pancreatische B-zellen Peptide: Kinetic and Konzentration von Proinsulin, Insulin and C-peptide in Plasma and Urine Probleme der Mezmethode Klinische und Literaturubersicht", *J Clin Chem Biochem*, **18**, 313-326 (1980).
- Boehm TM, Lebovitz HE, "Statistical analysis of Glucose and Insulin responses to intravenous tolbutamide; evaluation of hypoglycemic and hyperinsulinemic states", *Diabetes Care*, **479-490** (1979).
- National Committee for Clinical Laboratory Standards, "Procedures for the collection of diagnostic blood specimens by venipuncture: approved standards", 4th Ed, NCCLS Document H3-A4, Wayne PA, (1998).
- Turkington RW, Estkowski A, Link M, "Secretion of insulin or connecting peptide; a predictor of insulin dependence of obese diabetics", *Archives of Internal Med*, **142**, 1102-1105 (1982).
- Sacks BD, *Tietz Textbook of Clinical Chemistry*, 2nd Ed, Philadelphia, WB Saunders Co (1994).
- Kahn CR, Rosenthal AS, "Immunologic reactions to insulin, insulin allergy, insulin resistance and autoimmune insulin syndrome", *Diabetes Care*, **2**, 283-295 (1979).

Revision: 3 Date: 060712 DCO: 0653

Cat# 5825-300

Size	96(A)	192(B)
Reagent (fill)	A)	2ml set
	B)	1 (13ml) 2 (13ml)
	C)	1 plate 2 plates
	D)	1 (20ml) 1 (20ml)
	E)	1 (13ml) 2 (13ml)
	F)	1 (8ml) 2 (8ml)

For Orders and Inquiries, please contact



Tel: +1 949.951.2665 Email: info@monobind.com
Fax: +1 949.951.3539 Web: www.monobind.com

Please visit our website to learn more about our other interesting products and services.



EC REP CEpartner4U, Esdoornlaan 13, 3951DB Maarn, The Netherlands
www.cepartner4u.eu